

DATASHEET

AmbuMan® SAM is an inexpensive training manikin with all the well-known and sought-after Ambu look-listen-feel features.

The Ambu SAM incorporates the unique and patented hygienic system that enables training with a face piece and head bag.

CPR becomes the most important factor when a person suffers from cardiac arrest

Blood circulation of the body stops, including the perfusion of the brain, lungs and organs. Without treatment the patient will die. With chest compressions a blood circulation through the body is created.

Due to this, it is important that people are trained in this discipline.

Correct Compression and Open Airway

The AmbuMan SAM enables allows students to have a realistic ventilation experience. The compression indicator confirms the correct compression depth.

The closed torso shape with all the precise anatomical features allows students to fully learn the BLS technique, including the use of semi-automatic defibrillators (AEDs) and the positioning of the corresponding electrodes.

The unique hygienic system prevents cross contamination.

The correct compression depth and required ventilation volume can be monitored by the feedback indicator on the manikin.

AmbuMan® SAM



KEY BENEFITS

- Robust and reliable cost efficient BLS training manikin with anatomically correct closed Torso
- Visual feedback built into the manikin torso, provides instant feedback of compression amplitude
- Ambu hygienic system for mouth-to-mouth and mouth-to-nose ventilation
- Easy AED Training
- Compact manikin, can be easily carried around and stored because it is a lightweight manikin

Ambu

MATERIALS

Basic unit	ABS (Acrylnitril-Butadien-Styrol)
Head	Polyvinylchlorid, hart
Chest skin	Polyvinylchlorid, weich
Face piece	Polyvinylchlorid, weich
Headbag, disposable	Polyethylen
Trainingsuit	Polyester
Carrying bag	PVC-beschichtetes Nylon


SPECIFICATIONS

Weight	4,5 kg
Length	45 cm
Width	37 cm
Height	25 cm
Max. Ventilation volume	1800 ml
Max. compression depth	7 cm


ORDERING SPECIFICATIONS

Item no.


A268 003 000	AmbuMan SAM with 25 head bags and 1 face piece, supplied in a soft carrying bag
234 000 702	Head bags (sets of 100)
234 000 703	Face piece (sets of 5)

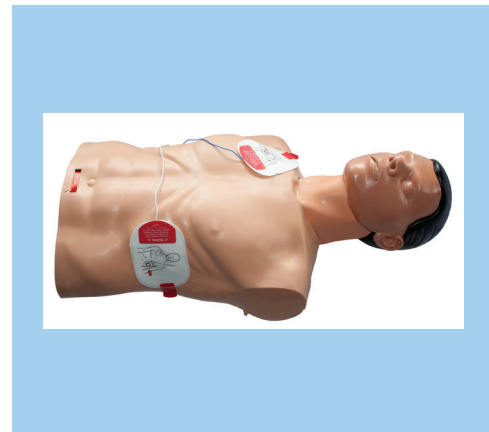

Placement of AED pads

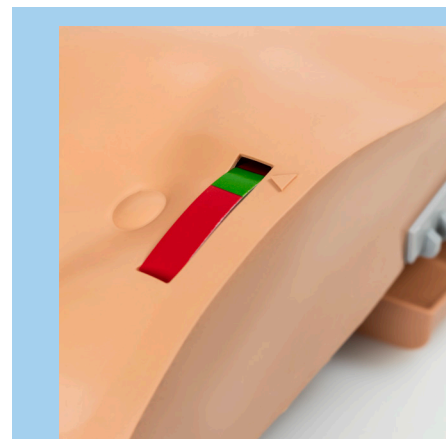

Compression depth indicator

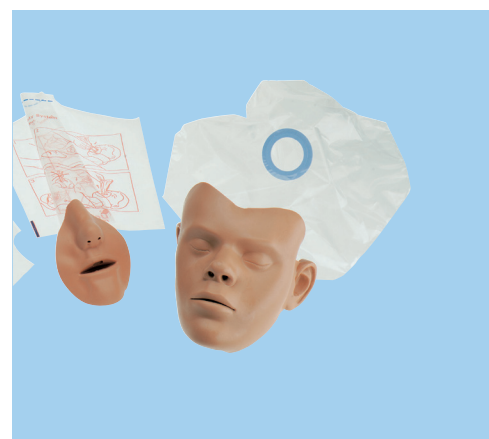

Hygienic System

Ambu

■ Ambu A/S
Baltorpbakken 13
2750 Ballerup
Denmark
T +45 72 25 20 00
F +45 72 25 20 50
ambu.com